



permobil



DYNAMIC BY DESIGN

Featuring a rear-wheel drive power base and two different speed packages of up to 6.5 MPH, the C350 Corpus 3G can easily keep pace with your busy lifestyle. Plus, the Corpus 3G seating system is designed for comfort beyond comparison, something you need in today's fast-paced world. The seating system follows the contours of the body, offers infinite adjustments for a custom fit and reclines 150°. The armrests feature built-in channels for cable routing, which greatly reduces the need for messy wire ties; while the legrests are curved for better support. Blending form and function, the C350 Corpus 3G is the ideal blueprint for mobility.

C350
CORPUS 3G



POWER FUNCTIONS



Power tilt and recline with elevating legrests

SHROUD COLORS



Crystal Blue




Vintage Wine



Platinum Silver



Xtreme Yellow

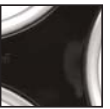

Black Onyx



50° of tilt with two range selections



150° power recline



8" power seat elevator



Power adjustable legrests with independently adjustable legrest height and footplate angle



4-wheel adjustable independent suspension



Independently adjustable armrest height, angle and pitch

SPECIFICATIONS

Maximum Speed	6.5 MPH	Seat Widths	17", 19", 21", 23" (Adjustable)
Driving Range	15-22 Miles	Seat Depths	14" - 22" (Adj. in 1" increments)
Max User Weight	300 lbs	Backrest Heights	20", 23" - 28"
Ground Clearance	3.1"	Seat To Floor Height	18.25"
Base Width	25"	Power Tilt	50°
Length w/o Front Riggings	36.5"	Power Recline	85° - 150°
Recommended Battery	Group 34	Mechanical Shear Reduction	5.25"
Gradient Capacity	10° Up 6° Dn	Power Legrest	85° - 170°
Power Seat Elevator	8"	Legrest Articulation Amount	8"
Arm Pad Options	10", 13", 16"	Electronics (standard)	R-Net
		Alternative Drive Controls	Yes w/R-Net



THE POWER OF MOBILITY

Permobil Inc.
300 Duke Drive
Lebanon, Tennessee 37090
800.736.0925
www.permobil.com